



Trulife  
Zumo Foot

 **steeper**  
Creating life's turning points, together


# Zumo Foot

The Zumo foot by Trulife is designed to withstand 205kg, making it ideal for users at the heavier weight level.

The split composite keel provides increased dynamic response, stability and improved medial/lateral compliance. The removable sandal toe foot shell comes in a choice of two colours, offers a high quality, natural cosmesis.

## Features and Benefits

- Heavy-duty foot withstands up to 205kg
- Durable composite design for up to K3 users
- Pre-assembled with a male pyramid for quick and easy installation
- Split composite keel for increased stability
- Dynamic toe and cushioned heel response
- Removable, life-like, sandal-toe foot shell

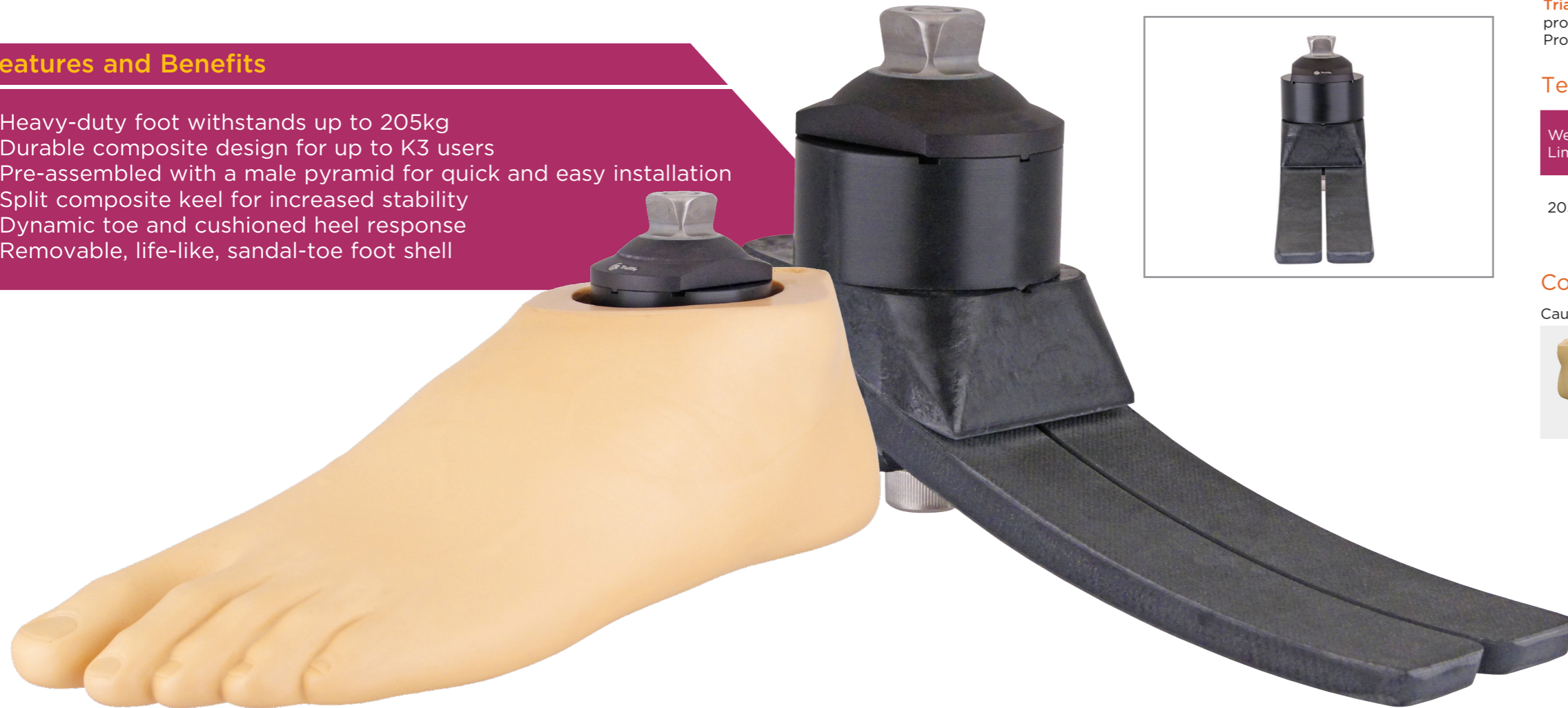

## Ordering Information

Part Number	Description	Available Sizes (cm)
NSZM450	Zumo Caucasian 9.5mm Heel Rise	22-30 cm
NSZM453	Zumo Dark Brown 9.5mm Heel Rise	22-30 cm

### Ordering Example - NSZM450-25-L

Part Number	Size	Left / Right
NSZM450	25	L

**Trial Period** - There is no trial period on this foot, however a test product may be available. If required please speak to your local Product Manager for more information.

## Technical Specification

Weight Limit	Impact Level	Weight of Unit	Sizes (cm)	Build Height (cm)	Warranty
205kg	L M H	750g Size 25	22-30	9 Size 25	2 years Foot shell - 6 months

## Cosmetics

Caucasian



Dark Brown



# 30mm Heavy Duty Structural and Socket Components

Trulife's range of heavy duty structural and socket components deliver both quality and dependability.

Made from industry graded titanium, these 30mm system components have a superior strength-to-weight ratio over the regular Aluminium components and a maximum weight limit of 205kg, making them ideal for use with the Zumo foot.

## Ordering Information

Heavy Duty Titanium Structural Components			
Part Number	Description	Weight Limit	
1	NAAAHD1210	4-Hole Male Pyramid Adapter	205kg
2	NAAAHD1217	Splined 4-Hole Male Rotatable Pyramid Adapter	205kg
3	NAAAHD1218	Splined 4-Hole Female Rotatable Receiver Adapter	205kg
4	NAAAHD1221	35mm Long Double-Ended Adapter	205kg
5	NAAAHD1225	Clamp Adapter	205kg
6	NAAAHD1200	Bonded Pylon 420mm	205kg
7	NAAASMP450	Socket Mounting Plate*	205kg

\*This includes fabrication components and Expulsion Valve as standard, but it can be used without, if a valve is not required.

### Weight Limit - 205kg

Please note that when ordering components from different ranges for the same build, the lowest rated component becomes the maximum weight limit of the build.

